

GlucoManager




Individual Management



Multiple Management

GlucoManager

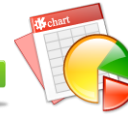

Patient Select

key in ID no.




View Reports

Download Data

Add/Edit

GlucoManager

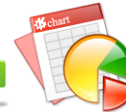

Patient Select

key in ID no.



Ann Lin
Abby Chan
Bob Lee
.
.
.
Zoe Wang



View Reports

Download Data

Add/Edit

GlucoManager



Profile

* must

Full Name* :

Advance

Gender : ☐ Male; ☐ Female

Date of Birth : YYYY-MM-DD, ☐ years old

ID No.* :

Medical Record No. :

Phone No. :

Address :

E-mail :

Meter Serial No. :

Previous

Reset

Save

Cancel

GlucoManager



Profile

Full Name* : Charlie Lo

Advance

Gender : ☐ Male; ☐ Female

Date of Birth : 1977-02-08, 35 years old

ID No.* : H122121666

Medical Record No. : 082

Phone No. : 03-668-6868

Address : 3F, No. 49-3, Lane 2, Sec. 2, Guang-Fu Rd.,
Hsin-Chu City, 30071, Taiwan

E-mail : charlie@epsbio.com.tw

Meter Serial No. : 27U003351

Previous

Reset

Save

Cancel

GlucoManager



Click “Advance” enter into following page

Schedule/Medications/ Regimen/Complications/Health Record/Inquiry

GlucoManager



Schedule

☒ Check this box to use different schedules for working and non-working days.

☐ AM/PM ☐ 24hr

☒ Meals ☐ Non-meal

Check boxes of non-working days:

Mon	Tue	Wed	Thu	Fri	Sat	Sun
<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input checked="" type="checkbox"/>	<input checked="" type="checkbox"/>

Working Day Schedule

Slot	From	To
Before-Breakfast:	上午 06:00	上午 09:00
After-Breakfast:	上午 09:00	上午 11:00
Before-Lunch:	上午 11:00	下午 02:00
After-Lunch:	下午 02:00	下午 05:00
Before-Dinner:	下午 05:00	下午 08:00
After-Dinner:	下午 08:00	下午 11:00
Night:	下午 11:00	上午 06:00

Non-Working Day Schedule

Slot	From	To
Before-Breakfast:	上午 06:00	上午 09:00
After-Breakfast:	上午 09:00	上午 11:00
Before-Lunch:	上午 11:00	下午 02:00
After-Lunch:	下午 02:00	下午 05:00
Before-Dinner:	下午 05:00	下午 08:00
After-Dinner:	下午 08:00	下午 11:00
Night:	下午 11:00	上午 06:00

Previous

Reset

Save

Cancel

GlucoManager



Medications

Insulin

- ☐ Glargine
- ☐ Lente
- ☐ NPH
- ☐ Premixed
- ☐ Rapid
- ☐ Regular
- ☐ Ultralente

- ☐ Other
- ☐ Custom 1
- ☐ Custom 2
- ☐ Custom 3

Oral Meds

- ☐ Acarbose
- ☐ Chlorpropamide
- ☐ Glibenclamide
- ☐ Gliclazide
- ☐ Glimepiride
- ☐ Glipizide
- ☐ Glipizide
- ☐ Glipizide
- ☐ Glipizide
- ☐ Glyburide
- ☐ Metformin
- ☐ Metformin/Glyburide
- ☐ Miglitol
- ☐ Nateglinide

- ☐ Other
- ☐ Custom 1
- ☐ Custom 2
- ☐ Custom 3

☐ Insulin Pump

☐ Diet & Exercise

Previous

Reset

Save

Cancel

GlucoManager



Regimen

Active/Regimen List

Regimen 1
Regimen 2
Regimen 3
Regimen 4
Regimen 5

Regimen Name: ☐ Active

Date Effective:

Meal Slot	Units	Type	Units	Type
Breakfast:	<input type="text"/>	<input type="text"/>	<input type="text"/>	<input type="text"/>
Lunch:	<input type="text"/>	<input type="text"/>	<input type="text"/>	<input type="text"/>
Dinner:	<input type="text"/>	<input type="text"/>	<input type="text"/>	<input type="text"/>
Night:	<input type="text"/>	<input type="text"/>	<input type="text"/>	<input type="text"/>


Previous

Reset

Save

Cancel

GlucoManager

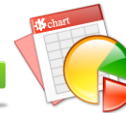

Complications

Complications:

- ☐ Cardiovascular
- ☐ Cerebrovascular
- ☐ Foot Infections
- ☐ High Blood Pressure
- ☐ High Cholesterol
- ☐ Nephropathy
- ☐ Neuropathy
- ☐ Peripheral Vascular
- ☐ Retinopathy

☐ Other

Previous

Reset

Save

Cancel

GlucoManager



Health Records

Diabetes Type	<input type="checkbox"/> Type I		
	<input type="checkbox"/> Type II		
Year Diagnosed	YYYY/MM/DD		
Height	<input type="checkbox"/> ???	<input type="radio"/> cm <input type="radio"/> ft.	
Weight	<input type="checkbox"/> ???	<input type="radio"/> kg <input type="radio"/> pound	
Before Meal	<input type="checkbox"/> --- <input type="checkbox"/>	<input type="radio"/> mg/dL <input type="radio"/> mmol/L	
After Meal	<input type="checkbox"/> --- <input type="checkbox"/>		
Hypo	<input type="checkbox"/>		
Ketones	<input type="checkbox"/> Negative	<input type="checkbox"/> Trace	<input type="checkbox"/> Small
	<input type="checkbox"/> Moderate	<input type="checkbox"/> Large	
Microalbumin	<input type="checkbox"/> Normal	<input type="checkbox"/> Positive	<input type="checkbox"/> >1000
	<input type="checkbox"/> ????		


Previous

Reset

Save

Cancel

GlucoManager

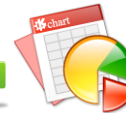

Health Records

Blood Pressure	Systolic	<u>???</u>	80 - 200
	Diastolic	<u>???</u>	40 - 150
Cholesterol	Total	<u>????</u>	⊙mg/dL ○mmol/L
	LDL	<u>???</u>	
	HDL	<u>???</u>	
Triglycerides		<u>????</u>	⊙mg/dL ○mmol/L
HbA1c		<u>??</u> %	4 - 15


Previous

Reset

Save

Cancel

GlucoManager

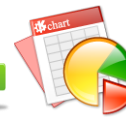

Inquiry

Date	Comment
Inquiry Note	
<input type="checkbox"/> Special View	


Previous

Reset

Save

Cancel

GlucoManager

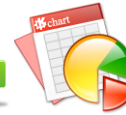

Patient : Charlie Lo (82)




Personal Data is saved!

OK

Cancel

GlucoManager

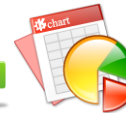

Patient : Charlie Lo (82)




Please confirm you are downloading data to
“full name (ID No.)

OK

Cancel

GlucoManager

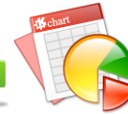

Patient : Charlie Lo (82)

Data is Downloading ...



Finish

Cancel

GlucoManager



Patient : Charlie Lo (082)

Unit of Measurement: ☒ mg/dL ☐ mmol/L

Logbook

Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Apply

Target: ■ low, ■ within, ■ high

Date ↓	Breakfast		Lunch		Dinner		Bedtime	
	Bef.	Aft.	Bef.	Aft.	Bef.	Aft.	Gluc.	Ins.
2012/03/12	89							
2012/03/13	92							
2012/03/14	87							
2012/03/14	90							
2012/03/15	100							
2012/03/16	93							
Average								
In Target								
# Results								

Scroll bar

Double click to "Record Details"

Record Details



Glucose

Date: 8/16/2006

Time: 11:56 AM

Slot: Before Lunch

Type: Finger Tip

Glucose: Or: ☐ HIGH

mg/dL

☒ Plasma

☐ Whole Blood

Glucose Comments

☐ Illness

☐ Food

☐ Exercise

☐ Exercise type

☐ Vacation

☐ Different Food

☐ Before

☐ Mild

☐ Stress

☐ Fasting

☐ During

☐ Hard

☐ Feel Hypo

☐ Too much food

☐ After

☐ Medication

☐ Menses

☐ Not enough food

☐ Other

Comments:

Reset

Save

Cancel

Record Details



Insulin

Date: 12/04/2002
Time: 04:20 AM
Slot: Night

☐ Bolus
 (Units)

☒ Injection
Type:
Amount: (Units)

☐ Pump Daily Total
 (Units)

... Add Insulin

Comments:

Reset

Save

Cancel

Record Details



Oral Meds

Date: 12/04/2002

Time: 04:22 AM

Slot: Night

Type:

Amount:

... Add Oral Meds

Comments:

Reset

Save

Cancel

Record Details



Exercise

Date: 12/04/2002

Time: 04:23 AM

Duration: hrs min

Intensity:

Comments:

Reset

Save

Cancel

Record Details



Food

Date:

Time:

Meal:

Carbs:

Fats:

Calories:

Proteins:

Comments:

Reset

Save

Cancel

GlucoManager



Patient : Charlie Lo (082)

Unit of Measurement: ☒ mg/dL ☐ mmol/L

Logbook

Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Apply

Target: ■ low, ■ within, ■ high

Date ↓	Breakfast		Lunch		Dinner		Bedtime	
	Bef.	Aft.	Bef.	Aft.	Bef.	Aft.	Gluc.	Ins.
2012/03/12	89							
2012/03/13	92							
2012/03/14	87							
2012/03/14	90							
2012/03/15	100							
2012/03/16	93							
Average								
In Target								
# Results								

Scroll bar

Double click to "Record Details"

GlucoManager



Patient : Charlie Lo (082)


Unit of Measurement: ☒ mg/dL ☐ mmol/L

Glucose Trend

Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Apply

Target: ■ low, ■ within, ■ high

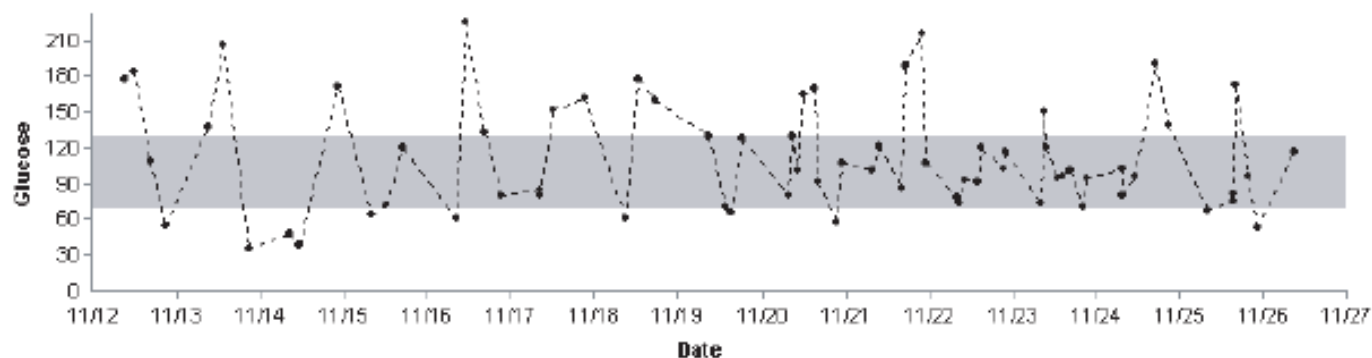

GlucoManager



Patient : Charlie Lo (082)

Unit of Measurement: ☒ mg/dL ☐ mmol/L


Pie Chart

Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Apply

Target: ■ low, ■ within, ■ high

Overall Totals

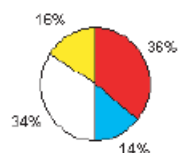


of Results

Above Target	25
Below Target	10
Within Target	24
Hypoglycemic	11

akfast


of Results

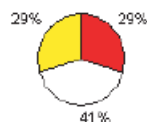

Above Target	5
Below Target	0
Within Target	7
Hypoglycemic	5

akfast


of Results

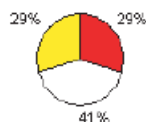

Above Target	5
Below Target	0
Within Target	7
Hypoglycemic	5

akfast


of Results

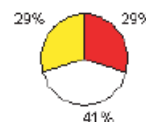

Above Target	5
Below Target	0
Within Target	7
Hypoglycemic	5

akfast

of Results

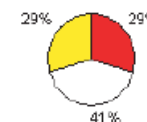

Above Target	5
Below Target	0
Within Target	7
Hypoglycemic	5

ikfast

of Results



Above Target	2
Below Target	1
Within Target	5
Hypoglycemic	0

ikfast

of Results



Above Target	2
Below Target	1
Within Target	5
Hypoglycemic	0

ikfast

of Results



Above Target	2
Below Target	1
Within Target	5
Hypoglycemic	0

ikfast

of Results



Above Target	2
Below Target	1
Within Target	5
Hypoglycemic	0

GlucoManager




Patient : Charlie Lo (082)

Standard Day

Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Target: ■low, ■within, ■high

Apply

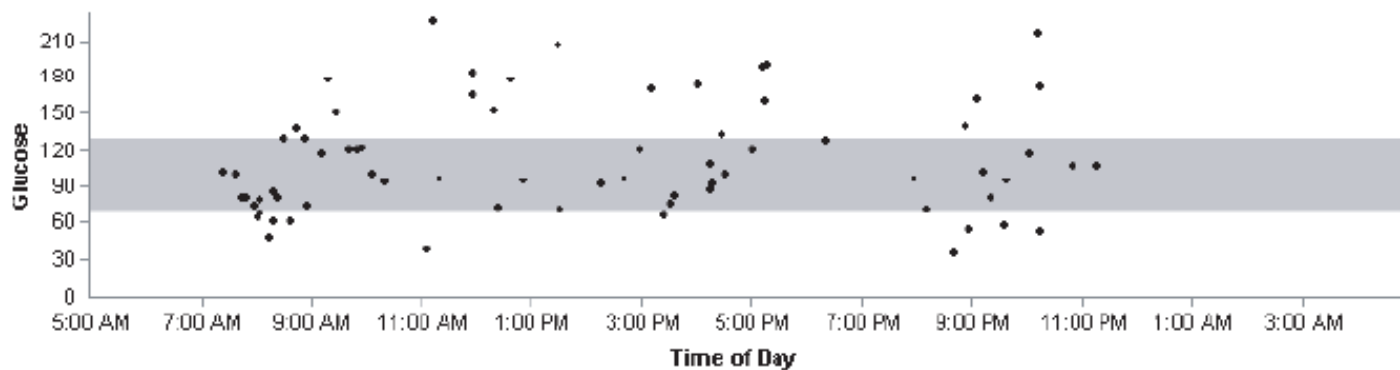

GlucoManager




Patient : Charlie Lo (082)

Average Readings


Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Target: ■ low, ■ within, ■ high

Apply

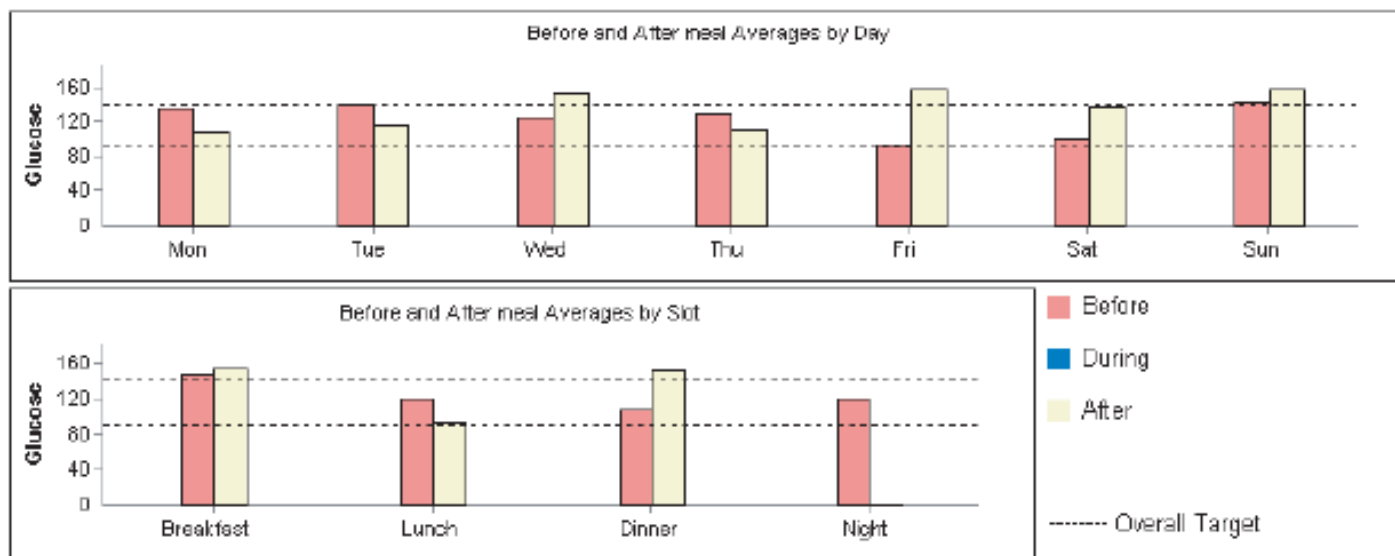

GlucoManager

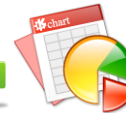

Patient : Charlie Lo (082)


Unit of Measurement: ☒ mg/dL ☐ mmol/L

Histogram

Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30


Apply

Target: ■ low, ■ within, ■ high

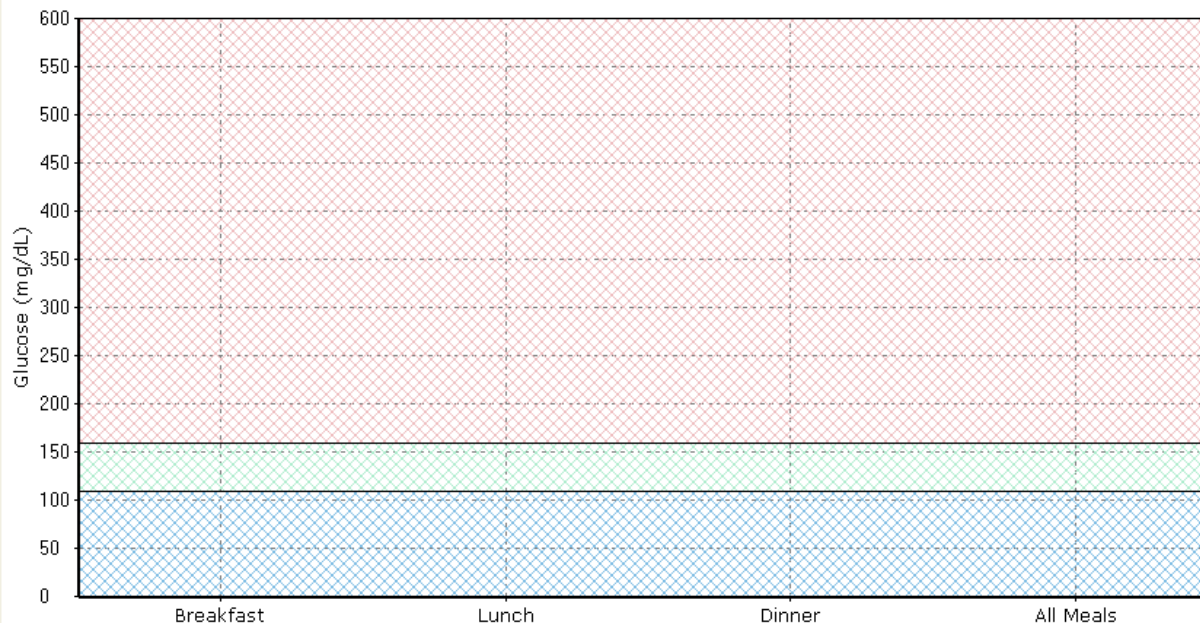

■ Pre-Meal average ■ Post-Meal average

GlucoManager

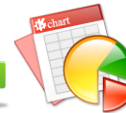

Patient : Charlie Lo (082)

Unit of Measurement: ☒ mg/dL ☐ mmol/L

Data List


Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Apply

Target: ■ low, ■ within, ■ high

Date	Time	Glucose	Insulin	Oral Med	Exercise	Food	Comments

Scroll bar

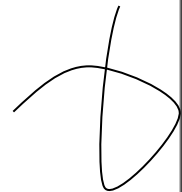

GlucoManager



Patient : Charlie Lo (082) Unit of Measurement: ☒ mg/dL ☐ mmol/L

Summary Checks

Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Apply

Target: ■ low, ■ within, ■ high

Full Name	: Charlie Lo				
ID No.	: H122121604				
Date of Birth	: 1977-02-08, 35 years old				
# Results	: 108			Standard Deviation: 8	
Above Target:	25	Below Target:	25	Within Target:	25
				Hypoglycemic:	32
Average	: 108				
Days Average	: 108 /7 days	: 108 /14 days	: 108 /30 days	: 108 /60 days	: 108 /90days
HbA1c	: 5.2%	Blood Pressure	: 110/79	Ketones	: Trace
Cholesterol	Total	: 197	LDL	: >500	HDL : >500

GlucoManager



Patient : Charlie Lo (082)

Unit of Measurement: ☒ mg/dL ☐ mmol/L

Summary Checks


Date Range: latest 7 days*
2012-03-01 ---- 2012-03-30

Apply

Target: ■ low, ■ within, ■ high

Full Name	: Charlie Lo
ID No.	: H122121604
Date of Birth	: 1977-02-08, 35 years old
# of results	: 108
Average result	: 108/7 days, 108/14 days, 108/30 days, 108/60 days, 108/90 days
Standard Deviation	:

GlucoManager

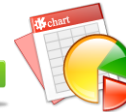


Multiple Patients Select

Download Data


Quick View

GlucoManager

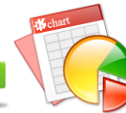

Medical Record No.	Patient Name	Glucose	Testing Date
082	Charlie Lo	89	2012-03-30
090	Abby Wang	210	2012-04-01
100	Nick Fong	115	2012-04-01
101		180	2012-04-03


Scroll bar

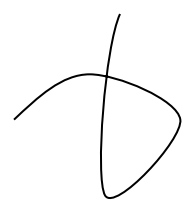

“Double click to Profile Page”


“Abnormal value will show color block”

GlucoManager

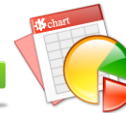

Patient	# of Readings	% Below Target	% Above Target	# Hypo	HbA1c	Pre-meal Average	Post-meal Average


Scroll bar

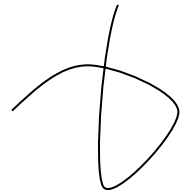

Select Patient:

[All Patients]
[Special View Patients]
Ann Lin
Bob Lee

GlucoManager

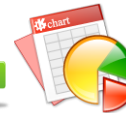

Archive

Retrieve

To File : C:\Program Files\LifeScan\OneTouch\Archive\Archive

Browse...

From : 10/14/2002

To: 11/14/2002

Select Patient :

[All Patients]


Doe, John

☐ Delete Readings After Archive

Archive

Cancel

GlucoManager

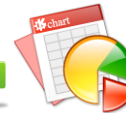

Archive

Retrieve

From File: C:\Program Files\LifeScan\OneTouch\Archive\Archive

Browse...

Examine File

From:


To:

Select Patient:

Archive

Cancel

GlucoManager

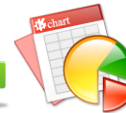

Export

Import

To File:

Select Patient: ▼

From: ▼

To: ▼

Archive

Cancel

GlucoManager



Export

Import

From File: C:\Program Files\LifeScan\OneTouch\Export\Export_

Browse...

Examine File

Select Patient:

From:

To:

Archive

Cancel

GlucoManager



Multiple Patients Select

Patient List

View Report

Library

GlucoManager



Patient : 082 (Charlie)



Are you sure to delete this Patient and his / her prescription?

OK

Cancel

GlucoManager



Are you sure to leave GlocuManager?

OK

Cancel